


FitRx™

Vibrating Foot Massager

User Guide

Welcome to the FitRx fam! Let's get you started using your new Vibrating Foot Massager:

- 1.** Place the foot massager on the floor.
- 2.** Position your foot on top of the massager.
 - a. Don't wear shoes while attempting to use the massager (socks are fine, though).
- 3.** Move your foot back and forth along the massager, controlling the pressure as you massage.
- 4.** If you wish to use the vibration function, simply press the power button on the side of the massager to turn it ON. Press the power button again to turn the vibration OFF.

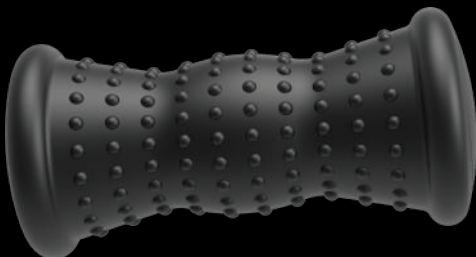

Features

- Textured surface for a deeper, more effective massage
- Massager can be used with or without vibration
- Ergonomic design to support natural foot arch

Benefits

- Soothing vibration stimulates blood flow
- Relieves muscle tension
- Decreases stress and promotes relaxation

relieve. recover. relax.

Disclaimer

Read the instructions carefully before using the massager. The use of this equipment is subject to following its instructions closely. tzumi is not liable for any damage or loss caused by inappropriate use of the device.

If you are: pregnant, diabetic, have a pacemaker, recovering from surgery, suffer from migraines and/or neuropathic disorders, disc protrusion and/or spinal inflammation, do not use this massager without the approval of your doctor. If you feel pain or discomfort while using the massager, stop using it immediately and consult your doctor.

Specifications

Size: 2.9 x 7.1 inch
Voltage: 3V

Speed: 4800R/min
Power: 2x AA batteries

Note: Open the plastic cover on the side of the massager to change the batteries.

**For more details on the
FitRx family of products
or to check out
more innovations visit
tzumi.com**



FitRX Vibrating Foot Massager, Distributed by tzumi Inc, NY, NY, 10016.
©Copyright 2021 tzumi Inc. All Rights Reserved. Made In China. All other trademarks and tradenames are the property of their respectful owners. All pictures shown are for illustration purpose only. Actual product may vary due to product enhancement. Subject to change without prior notice.