


VibraWrap

heated massage belt



User Guide

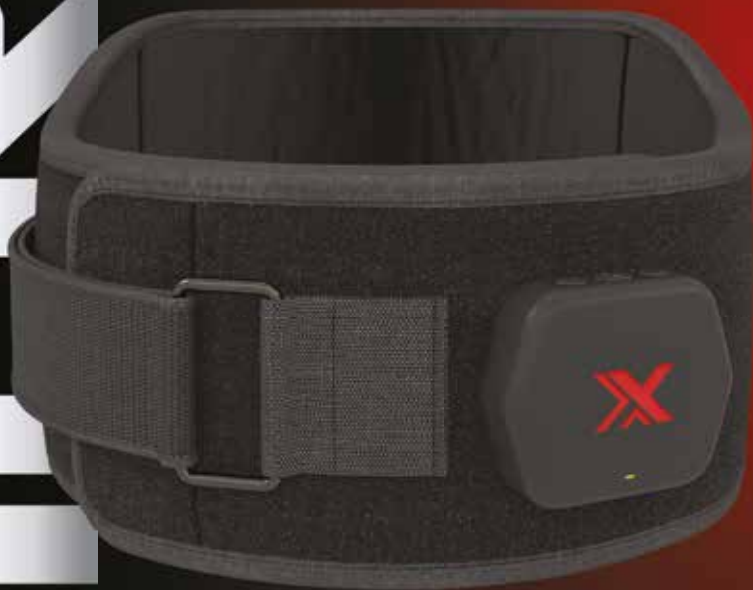

Congratulations!

You're now one step closer to improved physical fitness and overall wellness. The **FitRx® VibraWrap Belt** delivers results everyone can benefit from. Whether your goal is relief from workout fatigue or chronic joint pain, improved athletic performance and circulation, or a relaxing massage, vibration massage therapy is the perfect solution. We couldn't be happier to be taking this journey together while helping promote a healthier you from the inside out.

Welcome to the FitRX fam.

Are you ready?



Better Results For Your Fitness Goals.

What's Included

FitRx® VibraWrap Belt
Controller

USB-C Charging Cable
User Guide

Benefits

Increased Mobility & flexibility
Pre-workout warm-up
Post-workout recovery
Sore muscle relief
Improved blood circulation
Joint pain relief

Myofascial release
Trigger point massage
Relaxation

CUSTOMER SUPPORT

If you have any questions about setting up or using this product, we are here to help. Call or email us, and an experienced customer support representative will be happy to assist you: **1-855-60-TZUMI** or **support@tzumi.com**.



FitRxx

Total Body Fitness & Recovery

With the FitRx® VibraWrap, you can alternate between three heat levels (113°F, 122°F, and 131°F) to find the right temperature for you. Enhance your recovery experience with three vibration modes, providing up to 75 minutes of relief. With these customizable settings, you can reduce muscle fatigue, soreness, and improve blood circulation along your abdominal region.



User Instructions

Please read the following user instructions before using the FitRx® VibraWrap. Keep this user guide for future reference.

SETTING UP YOUR VIBRAWRAP

1. Before each use, always inspect the product to ensure there is no damage to the charging cable, the controller, or the vibrating belt.
2. Plug one end of the USB charging cable into the controller's charging port, and the other end into a USB power source. The LED indicator light will alternately flash red, blue and green while the unit is charging, and stop glowing once it is fully charged. Fully charge the VibraWrap prior to using it for the first time.
3. The LED indicator light will glow green when the unit is fully charged, glow blue when the battery level is at 50%, and glow red when it is at 30%. The LED indicator light will blink red when the unit is almost out of power.
4. Wrap the belt around your torso.

USING YOUR VIBRAWRAP

1. To turn the controller on, press and hold the Power button. The LED indicator light will glow white, indicating that the unit is in standby mode.
2. Press the Massage button to turn on the first massage mode. Press it again to cycle through the 3 massage modes.
3. Press the Heat button to turn the heat to level 1. Press it again to cycle through the 3 heating modes.
4. To turn the controller OFF, press and hold the Power button for 2 seconds. The LED indicator light will turn off.

NOTE: The VibraWrap Belt will turn off once 75 minutes have elapsed. To use the vibrating belt again, press and hold the Power button to turn the controller back on.

Cleaning & Care

1. Before cleaning the VibraWrap, ensure that the controller is turned off, and that it is not charging.
2. Use a damp, clean cloth to gently wipe both the controller and the belt.
3. Dry the unit thoroughly and store it in a cool, dry location, away from direct sunlight, when it is not in use.

Important Safeguards

If you have any of the following pre-existing conditions, please consult with your doctor before using the VibraWrap:

- Pacemaker or a similar implanted electronic device.
- Heart disease.
- Skeletal damage or fragility, especially in the spine or back.
- Recent surgery or operation.
- Currently or recently pregnant.
- Diagnosed with a malignant tumor.
- Osteoporosis.
- Issues with blood circulation or thrombosis.
- Skin infection of any kind.
- Abnormal blood pressure.

Warning: Consult with your doctor before using the VibraWrap on an injured part of the body.

Please stop using the VibraWrap immediately if any of the following occurs:

- The VibraWrap has been dropped or damaged in any way.
- The VibraWrap has been submerged in water.
- The VibraWrap is not operating correctly.

Safety Instructions

To reduce the risk of electric shock, fire, and injury, follow these instructions carefully.

- This device is recommended to be used by adults only.
- This device is for household use only.
- Frequent use of high pressure may cause abrasions. If you experience pain or discomfort, stop using the vibrating belt immediately.
- Give the vibrating belt time to cool down in between sessions.
- Do not use the vibrating belt in damp or humid environments, such as in the bathroom or near a pool.
- Exposure to water may cause the vibrating belt to short circuit.
- Exposure to high heat may cause the vibrating belt to overheat or explode.
- Do not use the vibrating belt during a thunderstorm.
- Do not cover the vibrating belt with anything while it is in use.
- Do not attempt to disassemble or modify the vibrating belt.
- Make sure the vibrating belt is turned off when setting up or cleaning it.



SPECIFICATIONS

Input: 5V/2A

Power: 45W

Battery Capacity: 7.4V/2800mAh

Charging Time: up to 4 hours

Working Time: up to 75 minutes

Product Size: 54.13 x 10.2 in.

Product Weight: 1.76 lbs. (0.8 kg)

Temperature Levels: 113°F, 122°F, 131°F

Disclaimer

Read this user guide carefully before using the massager. The use of this equipment is subject to following its instructions closely. Tzumi is not liable for any damages or losses caused by the inappropriate use of this device.

Do not use the massager without the approval of your doctor if you are: pregnant, diabetic, have a pacemaker, recovering from surgery, suffer from migraines or neuropathy disorders, disc protrusion or spinal inflammation, do not use this massager without the approval of your doctor. If you feel pain or discomfort while using the massager, stop using it immediately and consult your doctor.

FitRx® VbraWrap, item No. 9487, Distributed by tzumi Inc. NY, NY, 10016. ©Copyright 2024 tzumi Inc. All Rights Reserved. Made in China. All other trademarks and trade names are the property of their respective owners. All pictures shown are for illustration purposes only. Actual product may vary due to product enhancement. Subject to change without prior notice.

For more products and innovation
visit FitRxRecovery.com.

FitRx